

**REMARKS**

Claims 1-22 are currently pending in the application. Claims 1-22 were rejected. Claim 1 has been amended.

The Examiner objected to the Information Disclosure Statement filed on June 3, 2004, as not complying with the requirements of Rule 98. The requested information was resubmitted in a communication dated December 16, 2004, and the objection is believed addressed thereby.

The Applicant respectfully requests full consideration of the art submitted.

The Examiner indicated that the claim to prior applications on the first page of the application should be updated to include reference to the recently issued patent corresponding to application to which that section refers. The Examiner also identified a typographical error on page 22 of the application. The specification has been amended as requested and the objections are believed addressed thereby.

The Examiner rejected claims 1-22 under 35 U.S.C. 112, first paragraph, for failing to comply with the written description requirement. Specifically, the Examiner referred to the limitation in claim 1 which recites "the form factor corresponding to a thickness of the cassette enclosure which is less than about 15 mm," as not being described in the specification. The Examiner indicated that, to the extent such a limitation is described in IEC 60406 (incorporated by reference), the specification should be amended to include such description. The rejection is respectfully traversed.

The Applicant directs the Examiner's attention to Figs. 9 and 22-24, each of which shows a cassette enclosure having a thickness of either 0.59" or 0.6". As is well known, this thickness corresponds to about 15 mm, i.e., 25.4 mm to an inch. The Applicant therefore respectfully submits that the disclosure as filed supports the claim limitation to which the Examiner objects.

Notwithstanding the foregoing, the specification has been amended as detailed above to make specific reference to the subject matter to which the Examiner referred. Specifically, the

specification has been amended to augment these explicit references in the drawings with information from the third edition of IEC 60406 (in force on the priority date of the present application) regarding the dimensions of standard size radiography cassettes (see page 13 of that document). As requested by the Examiner, the undersigned hereby declares that no new matter is introduced by this amendment, in view of which, the rejection of claims 1-22 under section 112 should be withdrawn.

It should also be noted that claim 1 has been amended to more clearly correspond to the requirements of IEC 60406. That is, claim 1 now recites that the cassette thickness has “a maximum of about 15 mm.” No new matter is introduced by this amendment. In addition, this amendment does not require further search in that it is not significantly different from the previous limitation in this regard.

The Examiner rejected claims 1, 9, 11, 12 and 20 under 35 U.S.C. 103(a) as being unpatentable over International Publication No. WO/99/28765 (Mueller) in view of U.S. Patent No. 5,221,843 (Alvarez) and the third edition of IEC 60406. The Examiner also rejected claim 10 as being unpatentable over Mueller, Alvarez, and IEC 60406, and further in view of U.S. Patents No. 6,239,516 (Floresta) and No. 5,912,944 (Budinski). The Examiner also rejected claims 13-15 as being unpatentable over Mueller, Alvarez, and IEC 60406, and further in view of U.S. Patent No. 5,757,021 (Dewaele). The Examiner also rejected claims 16-19 as being unpatentable over Mueller, Alvarez, and IEC 60406, and further in view of U.S. Patent No. 5,864,146 (Karellas). The Examiner also rejected claim 21 over Mueller, Alvarez, and IEC 60406, in view of Floresta, and claim 22 over Mueller, Alvarez, and IEC 60406, in view of Applicant’s admitted prior art. The rejections are respectfully traversed.

Mueller describes an x-ray “cassette” which contains a phosphor carrier and a device for reading out information stored in the phosphor carrier. The device includes a radiation source for exciting the phosphor carrier, and a receiving means for receiving the resulting radiation. See

Abstract. Referring to the description in U.S. Patent No. 6,373,074 (i.e., the U.S. patent claiming priority from Mueller), the receiving means is described as a CCD 12 in combination with a Selfoc lens assembly 14 which focuses the excitation radiation on the individual elements of the CCD. See column 5, lines 10-27. Notably, Mueller states that “[i]t is possible to limit the thickness of the x-ray cassette to about 45 mm.”

As an initial matter, it is important to note that the term “cassette” as used by Mueller should not be equated with a “standard radiographic film cassette” as recited in the claims of the present application. That is, as would be apparent to those of skill in the art, the “cassette” to which Mueller refers is more commonly referred to as a “bucky,” which is a term used in the industry for the cassette tray which also includes a reciprocating grid above it (invented by Bucky). Like radiographic film cassettes, buckys also have standard sizes (although not subject to an international standard) and the typical clearance for a cassette tray is 1.94” (about 49 mm). Background information relating to the nature of buckys has been provided herewith to illustrate this important distinction. This information includes illustrations created by the Applicant to facilitate the Examiner’s understanding of the difference between a “bucky” and a cassette as commonly understood in the art.

Despite the language in Mueller and Alvarez regarding compatibility with conventional equipment, neither of the references teaches a device which can conform to the form factor recited in claim 1 and referred to in the present specification. Rather, Mueller indicates that its lower limit on “cassette” (i.e., bucky) thickness is at least three times the upper limit of the preferred standard cassettes defined by IEC 60406. This is due to the fact that the receiving means in Mueller, i.e., the CCD and Selfoc lens system, cannot be compressed below the stated limit due in large part to optical considerations. This is to be contrasted with the present invention which offers a variety of embodiments which are not so constrained.

Therefore, because Mueller’s system cannot be enclosed in a cassette having a “form

factor corresponding to a thickness of the cassette enclosure having a maximum of about 15 mm" as recited in claim 1, the rejection of claim 1 over the combination of Mueller, Alvarez, and IEC 60406 is believed overcome. Indeed, the fact that the apparatus taught by Mueller cannot fit within the maximum thickness prescribed by IEC 60406 makes it clear that the combination of Mueller's teachings with those of IEC 60406 is improper. The rejections of claims 2-22 are believed overcome for at least the reasons discussed.

The Examiner disagreed with the Applicant's arguments in the previous response stating that the 45 mm limit to which Mueller refers "is merely an example of the very small dimensions of the x-ray cassette and is not an express teaching of a lower limit of ~45 mm." The Examiner went on to say that the "applicant's argument that the cassette thickness cannot be manufactured with very small dimensions due to optical considerations such as the Selfoc lens is not persuasive since the Selfoc lens is optional."

With regard to the first point, the Applicant respectfully disagrees. In describing a key advantage, Mueller states that his technique makes it possible "to limit the thickness of the x-ray cassette to about 45 mm." It strains credibility to assert that assert that this is "not an express teaching of a lower limit." In addition to use of the term "limit" with reference to the phrase "about 45 mm," the inventor is describing a key advantage of his invention. To assert that the inventor would not refer to the absolute minimum value he thought possible at the time is simply not an accurate reading of the reference.

With regard to the second point, the Examiner mischaracterizes Mueller's statement at column 5, lines 12-14. Mueller is not saying that some form of lens is not required, but that having an individual Selfoc lens "for each stimulable point of the line of the phosphor plate...is not required for the invention." That is, the large array of Selfoc lenses contemplated by lines 12-13 of column 5 is not necessary. However, as would be understood by one of ordinary skill in the art, Mueller does not suggest that optical elements (i.e., a Selfoc lens or the equivalent)

between the plate and the photodetector are unnecessary. To the contrary, it would be apparent to one of ordinary skill in the art that Mueller's system would not work without some form of relay optics to focus the radiation from the plate on the CCD.

That is, the manner in which Mueller's technique reads out the information stored in a phosphor plate and the size of the CCD requires some distance between the phosphor plate and the CCD assembly which, in turn, requires an optical system to transmit light from the plate to the CCD. This can be understood with reference to Fig. 1 of Mueller.

As shown in that figure, the size of CCD assembly 12 and the requirement that it be placed at an angle and off to the side of laser diodes 11 (so as to be able to receive the emitted radiation) requires a minimum physical spacing which, in turn, requires the use of an intervening optical relay. This, therefore, forces a minimum size on the assembly which Mueller himself states to be "about 45 mm." The Applicant has provided some background information herewith (Schaetzing et al.) which further illustrates the necessity of intervening optics in a design like Mueller's.

By contrast, the readout mechanism employed by the present invention allows the image capture devices (e.g., photodetector array 114 of Fig. 1) to be placed in extremely close proximity to the surface of phosphor plate 104. This allows for embodiments without bulky intervening optics which, in turn, enables a much thinner assembly geometry (e.g., 15 mm or less) to be achieved.

Thus, because of the geometry imposed by Mueller's readout technique, a device constructed in the manner described cannot even approach a "form factor corresponding to a thickness of the cassette enclosure having a maximum of about 15 mm."

In view of the foregoing, the teaching of Mueller cannot be combined with IEC 60406 to obviate any of the claims of the present invention. The Applicant therefore respectfully requests that the rejection of claim 1 over Mueller, Alvarez, and IEC 60406 be withdrawn. In addition,

dependent claims 2-22 are also believed to be allowable over the cited art for at least the reasons discussed.

In view of the foregoing, Applicants believe all claims now pending in this application are in condition for allowance. The issuance of a formal Notice of Allowance at an early date is respectfully requested. If the Examiner believes a telephone conference would expedite prosecution of this application, please telephone the undersigned at (510) 663-1100.

Respectfully submitted,  
BEYER WEAVER & THOMAS, LLP



Joseph M. Villeneuve  
Reg. No. 37,460

P.O. Box 70250  
Oakland, CA 94612-0250  
(510) 663-1100




**Materials prepared by Michel Sayag  
in support of response to office action  
dated January 11, 2005**




### Size difference between Muller and Sayag designs




Muller design (per Agfa publication)

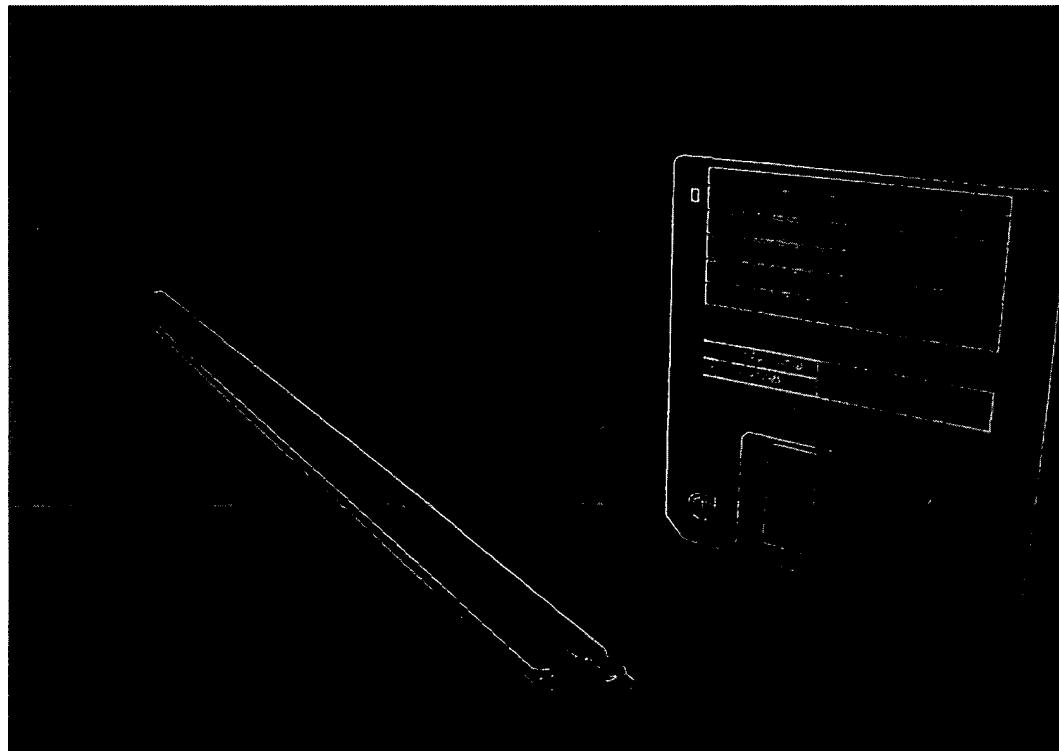

Sayag design

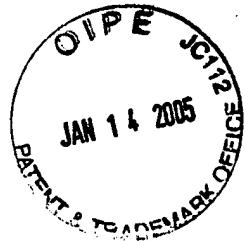

## X-ray Machine Components

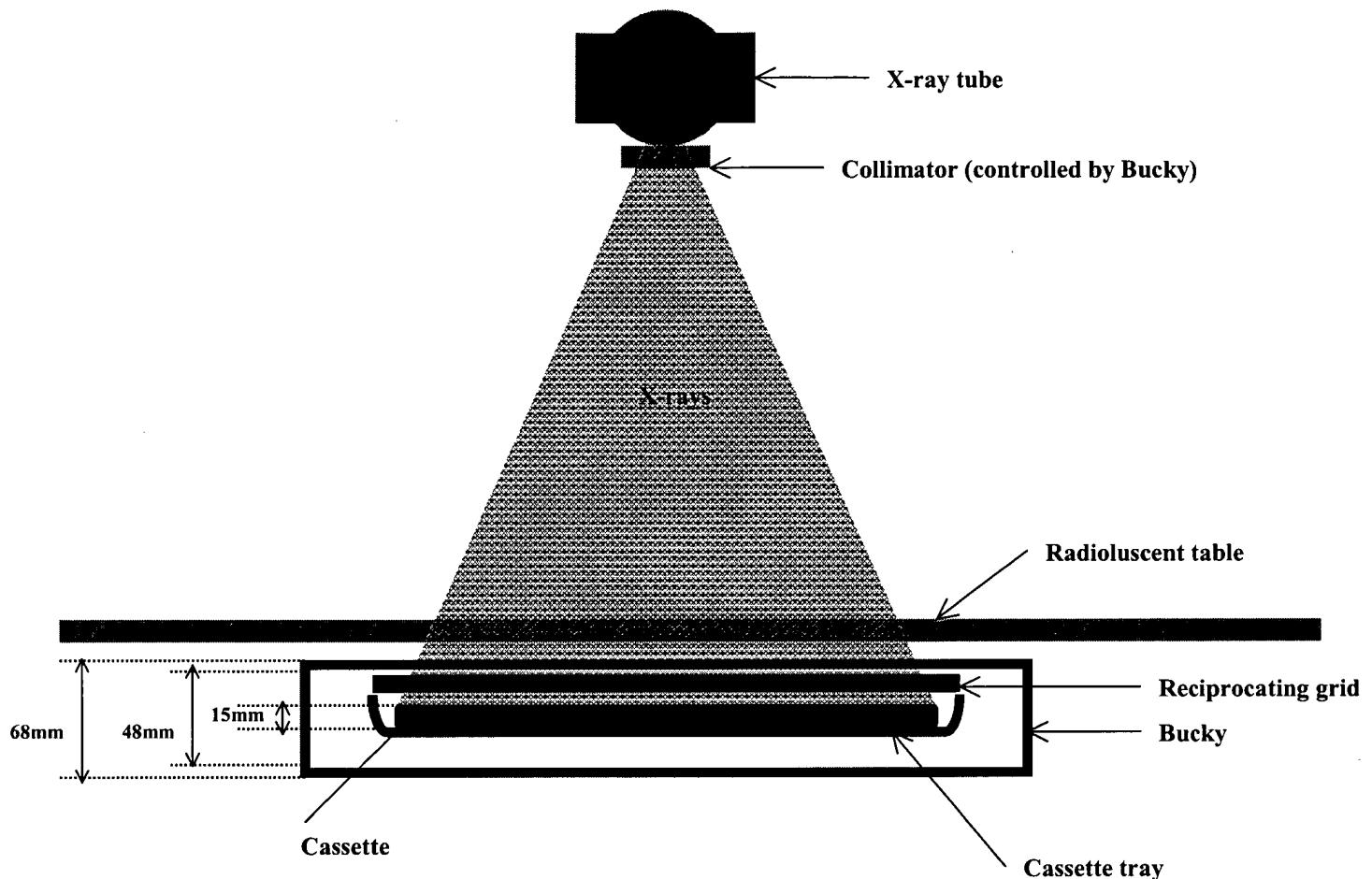


Bucky dimensions per Poersch drawing ([http://www.poerschmetal.com/Bucky\\_folder/Bucky.html](http://www.poerschmetal.com/Bucky_folder/Bucky.html))



## Digital Bucky versus Digital X-ray Cassette

Muller enclosure (Bucky) insertable in existing x-ray units

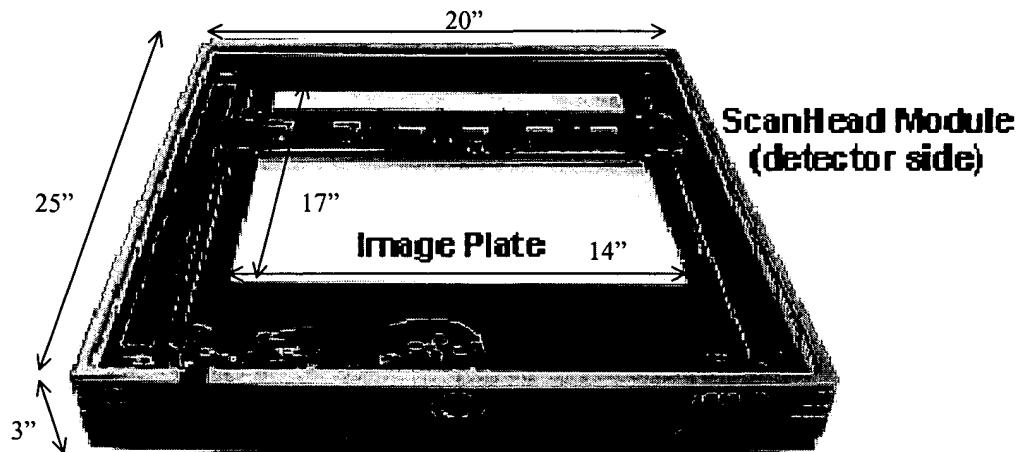


**Figure 12.** View inside a CR flat panel containing the ScanHead® module. The stimulation section is underneath the image plate.

Agfa's original text relating to Fig.12:

"Fig. 12 shows a prototype of a ScanHead®-based flat-panel CR system that might be used for portable radiography or as the capture module inside a Bucky table."

Sayag enclosure (cassette) insertable in cassette tray of existing x-ray units

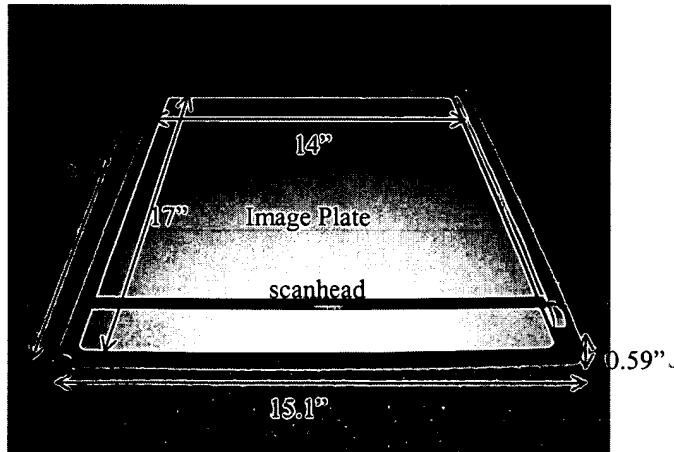

Dimensions taken from Sayag's Fig. 23 (no new information added)

# New high-speed scanning technique for Computed Radiography

Ralph Schaetzling<sup>1</sup>, Robert Fasbender, Peter Kersten  
Agfa-Gevaert AG., HealthCare, Munich, Germany

## ABSTRACT

The natural luminescence decay time of storage phosphors limits the scan speed of today's flying-spot Computed Radiography (CR) scanners, since this "afterglow" emission degrades the MTF along the fast-scan direction. Although this sharpness loss can be compensated to some extent with image processing, there is a decrease of signal-to-noise in the process. Higher scan speeds can also decrease the total amount of signal captured at each point, that is, the read-out depth. One may reduce these problems by scanning entire lines at once, rather than points. However, this can create other problems related to stimulation and collection efficiency, system gain, sharpness, and noise. We have developed a new scanning engine containing a linear, laser-diode-based stimulation source and a light collection/detection system with specially designed optics and multiple, linear, asymmetric CCDs. The stimulation/collection systems are housed in a compact, movable read-head that scans quickly (43x43 cm in 5 s, independent of resolution) over a stationary image plate (IP). We evaluated this so-called ScanHead® system with conventional powder IPs and with new CsBr:Eu<sup>2+</sup> needle-based IPs. We also compared it to a conventional flying-spot scanner. Image quality results differed depending on plate type. DQE values for conventional powder IPs in the prototype scanner were comparable to today's state-of-the-art commercial CR systems (DQE(0)~0.25-0.28). DQE values for prototype needle IPs were significantly better than for powder IPs (DQE(0) ~0.58-0.60) and comparable to those quoted for flat-panel DR systems based on CsI:Tl.

**Keywords:** Computed Radiography, CR, storage phosphor, digital radiography, DR, image quality, scanning, needle image plate, IP, CsBr

## 1. INTRODUCTION

Although some of the technologies needed to do Computed Radiography (CR) have existed for over 60 years<sup>1</sup>, the first commercial CR system was introduced only about 20 years ago<sup>2</sup>. Today, CR is a clinically accepted imaging modality for practically every exam type, and many manufacturers have entered the market with CR systems and related equipment. Over the last 20 years, manufacturers have been able to decrease the physical size and increase the throughput of CR systems considerably, while increasing their image quality. However, there are some inherent limitations to these tendencies.

All commercial CR scanners today use, more or less, the principle of flying-spot scanning, in which a tightly focused laser beam is used to stimulate the latent image in a moving storage phosphor screen (image plate) **one point at a time** over the entire screen surface (Fig. 1). Through appropriate collection optics, some (preferably large) fraction of the emitted light from each point is captured, converted by a photodetector (usually a photomultiplier tube) into an analog electrical signal, then sampled and quantized to produce a digital image. The hardware needed to produce and deflect the laser beam, as well as the components needed to collect the emitted light and convert it into an electrical signal require a certain amount of space, which makes it difficult to reduce indefinitely the physical size of flying-spot scanners. In addition, these discrete components add cost and complexity.

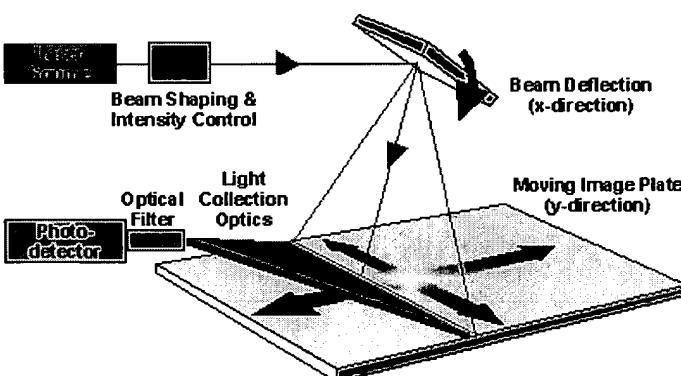


Figure1. Schematic diagram of the components in a flying-spot scanner.

<sup>1</sup> [ralph.schaetzling.rs@germany.agfa.com](mailto:ralph.schaetzling.rs@germany.agfa.com); Phone:+49(89)6207-3352; FAX:+49(89)6207-7619;  
<http://www.agfa.com/healthcare>; Agfa-Gevaert AG, Tegernseer Landstr. 161, 81539 Munich, Germany

The throughput of a flying-spot CR scanner depends on many factors<sup>1</sup>, one of them being how much of the stored latent image in the IP the scanner must/can extract in order to make an output image with reasonable image quality (this signal, a component of *system gain*, largely determines SNR). The fraction of stored signal released from the IP during scanning, called the *read-out depth*, depends on the amount of stimulation energy deposited per unit screen area, which is proportional to the incident laser intensity and exposure time. This is illustrated in Fig. 2<sup>3</sup>. Increasing the exposure increases the total emitted signal until the stored latent image is essentially depleted (any remaining signal must be erased before the IP can be reused). The read-out depth per pixel depends on the amount of time the laser beam spends "over" the pixel area, also known as the *dwell time*. Dwell time is inversely proportional to beam velocity (see Fig. 3).

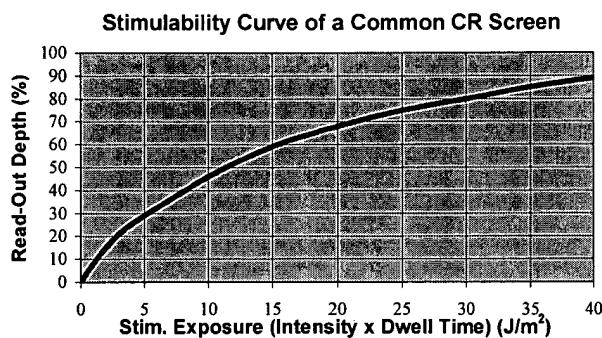


Figure 2. Stimulability curve for a common storage phosphor screen<sup>3</sup>. The fraction of stored signal read out increases with stimulating exposure until the latent image has been essentially depleted (any remaining signal on the IP must be erased before the next use).

So, the fundamental throughput trade-off in a flying-spot scanner is scan time vs. signal. Higher scan speeds reduce scan time and, therefore, increase throughput, all other things remaining constant. However, higher scan speeds also reduce dwell time, which lowers read-out depth and decreases image quality. From Fig. 3, one can see that increasing the laser intensity can compensate for this exposure reduction at higher velocities. However, increasing the laser intensity also broadens the incident beam within the phosphor layer (due to light scattering) which reduces MTF and produces a loss of image quality.

An additional throughput trade-off with flying spot scanners is "afterglow." If the stimulating laser beam exposing a storage phosphor screen were turned off suddenly, the screen would continue to emit light for a short, material-dependent period of time, called the *luminescence decay time*,  $\tau$ . The luminescence decay times of a few common storage phosphor materials,  $\text{BaFBr}_{x}\text{I}_{1-x}\text{Eu}^{2+}$  being the most commonly used today, are listed in Table 1<sup>4,5,6</sup>. If the dwell time gets short enough that it begins to approach this natural decay time, the light being collected at the current beam position contains "leftover" light still coming from the previous beam position(s). As a result, the image signal is effectively smeared in the fast-scan direction, lowering the MTF in that direction. This puts an upper limit on the throughput of flying-spot scanners (modern CR flying-spot scanners have dwell times on the order of 1-6  $\mu\text{s}$ , i.e., about 1.5-8 $\tau$  for a  $\text{BaFBrI}$  screen). Since the signal decay/blur characteristics are usually known, one may use image restoration techniques<sup>7</sup> effectively to undo this smearing effect, but only at the expense of lower SNR in the final image.

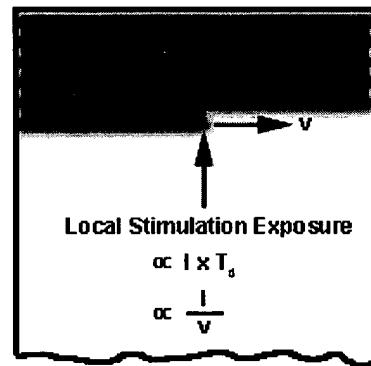


Figure 3. Schematic diagram of an IP being scanned in a flying-spot scanner. The local exposure is proportional to the intensity,  $I$ , and the dwell time,  $T_d$ , of the laser in that area.  $T_d$  is inversely proportional to the beam velocity,  $v_f$ .

Storage Phosphor	Luminescence Decay Time, $\tau$
SrS:Ce,Sm	~25 ns
RbBr:Tl	350 ns
CsBr:Eu <sup>2+</sup>	700 ns
$\text{BaFBr}_{x}\text{I}_{1-x}\text{Eu}^{2+}$	~750 ns

Table 1. Luminescence decay times (at room temperature) for several storage phosphor materials<sup>4,5,6</sup> ( $\text{BaFBr}_{x}\text{I}_{1-x}\text{Eu}^{2+}$  is the most commonly used in today's CR systems)


As illustrated above, the image quality of a CR system is a complex combination of, and often a trade-off between, material properties, hardware design, and processing software. Nevertheless, the image quality of CR systems has improved considerably over the last twenty years. Advances in storage phosphor chemistry, screen design, and manufacturing techniques have contributed especially to this improvement. If we use the increasingly accepted criterion of Detective Quantum Efficiency<sup>8</sup> (DQE) as a measure of image quality, then this improvement can be made more quantitative.

The DQE results reported for first-generation, flying-spot CR systems<sup>9</sup> were modest compared to existing screen/film systems. The DQE value at an exposure of 1 mR ( $10\mu\text{Gy}$ ) and a spatial frequency of 0 cy/mm, DQE(0), was about 11%. At 2 cy/mm, that value went down to around 5%. This tended to limit the clinical uses of these initial CR systems to applications that benefited from CR's wide exposure latitude, but that were less dose- and resolution-sensitive (e.g., portable radiography). An interesting study<sup>10</sup> of four more recent generations of CR systems (scanners and IPs) reported that the cumulative, relatively small improvements in scanner hardware and IPs had increased DQE(0) values by approximately 30%, while DQE values at higher frequencies had almost tripled. Thus, CR systems had become clinically competitive with screen/film systems for many more diagnostic applications, even as screen/film systems had improved their DQE performance over time. Current DQE(0) values for state-of-the-art flying-spot CR systems are in the range of 25-28%. A more recent study<sup>11</sup> claims DQE(0) values that are 30-40% higher than current values through the use of a transparent support material for the IP. By reading the phosphor layer from both sides simultaneously (i.e., dual collection optics and dual photodetectors), both the signal and SNR are improved (the MTF does not change significantly, since the light scattering in the screen remains and the signal transmitted through the screen has little high-frequency content).


## 2. THE NEW SCANNING ENGINE

A number of the above limitations shared by flying-spot scanners can be overcome by using a novel scanning technique in which the IP is scanned **one line at a time**, rather than one point at a time. In particular, the size limitations imposed by the bulkier flying-spot components are no longer relevant. Furthermore, the trade-offs between throughput and image quality can be largely circumvented. For example, dwell times in the new technique are measured in milliseconds rather than microseconds, yet scan times are still significantly shorter than those of flying-spot scanners. In fact, with the new needle-based CsBr:Eu<sup>2+</sup> storage phosphors, the new scan engine produces image quality (DQE) that is comparable to that of the latest flat-panel DR systems based on CsI:Tl and amorphous Si flat-panel arrays<sup>3, 12</sup>.

Figure 4 shows a diagram of the concept and Figure 5 shows an actual prototype of the new scanning engine, dubbed ScanHead®. The laser sources, beam-shaping optics, intensity control, light collector, filter, and photodetectors are all contained in one compact unit that scans quickly (43 x43 cm in approx. 5s) over a stationary IP. The control electronics to pre-process the captured signal are also on board, so that the raw, 14-bit digital image data can be transferred directly from the ScanHead® to a control computer for further processing and distribution. In the following sections, we will describe some of the important design considerations for several of the major components of this new scan engine.

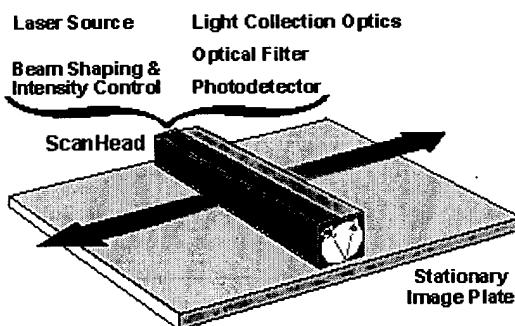

**Figure 4.** Schematic diagram of the new scanning concept. The compact, movable ScanHead® module contains lasers, beam shaping optics, intensity controllers, a light collector, an optical filter, and a photodetector.

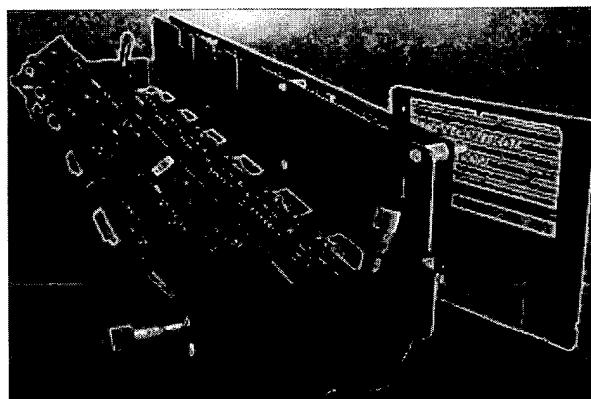


**Figure 5.** Laboratory prototype of the ScanHead® module.


## 2.1 Laser source, beam shaping, and intensity control

The light source is used to stimulate the storage phosphor locally to emit its stored energy. Although much of the installed base of CR systems still uses gas lasers for this task, newer CR systems are all based on laser diodes. These are available at wavelengths (red) and output powers (tens of mW) that are well matched to the stimulation spectra and sensitivities of today's storage phosphor materials. In addition, they are very compact light sources.

For ScanHead®, where an entire line must be illuminated at once, a linear array of laser diodes is used (Fig. 6). Although the sampling resolution of the CCD (Charge-Coupled Device) photodetector is 50 µm (see Sect. 2.4), providing one laser diode for each detector element would be cost-prohibitive and unreliable. Instead, we use the divergent, strongly asymmetric output beam profile of laser diodes to reduce the number of emitters to several dozen. The diodes are placed so that the strongly diverging (horizontal) beams overlap along the ScanHead® axis, producing a continuous line of stimulating light with sufficient power to reach an acceptable read-out depth in the IP. By inserting two plano-convex cylindrical lenses into the beam path, it is possible to focus the less diverging (vertical) beam axis down to 80 µm transverse to the line to maintain a reasonable MTF in that direction (as already indicated, the actual line width within the screen is largely determined by light scattering).

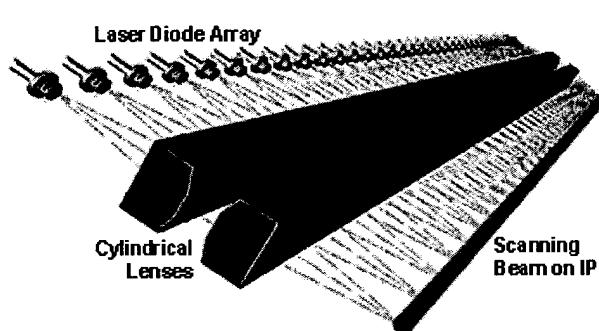

**Figure 6.** Schematic diagram of the linear laser diode array and the optics used to stimulate each line on the image plate. The large beam divergence along the lens axis creates the continuous scanning line. The smaller beam divergence across the lens (the rays in the figure) is brought to focus at the IP surface.

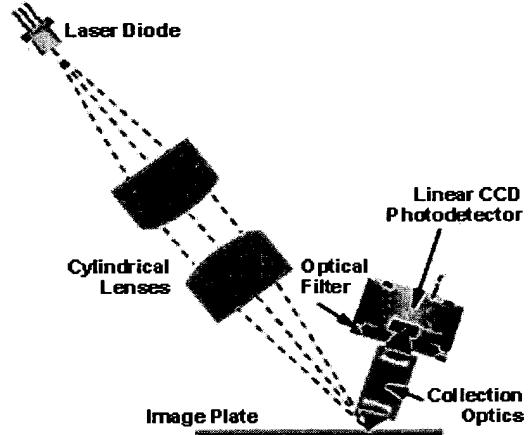

**Figure 7.** Cross-section of ScanHead® showing the angulation of the stimulation and detection sections.

Intensity control is needed to minimize static and dynamic exposure fluctuations to the IP. Such fluctuations easily can cause visible banding artifacts in the image. Static intensity variations along the scan line can occur, for example, from variability in the laser diode outputs and from overlap of the beam profiles of neighboring diodes. As in flying-spot scanners, such static intensity variations can be calibrated out by measuring the scanned line profile for a flat-field exposure and correcting all subsequent scanned data according to the measured profile. The tolerance for variability is fairly large here, on the order of 10-15%. The tolerance for dynamic fluctuations is considerably tighter. With an operating point on the linear part of the stimulability curve (Fig. 2), the allowable intensity variations may be less than fractions of 1%. ScanHead® operates near the saturated part of the stimulability curve (high read-out depth), which lowers this tolerance to around 1%. Correction of such dynamic intensity fluctuations must be done using active feedback, which can be done fairly straightforwardly in laser diodes, since they are driven directly.

Because of the compact nature of the ScanHead® module, it is difficult to achieve the usual beam geometry of a flying-spot scanner (stimulating beam nearly perpendicular to the IP). For the reflective configuration we have assumed so far, the stimulation and detection sections in ScanHead® are on the same side of the IP and must, therefore, be angled relative to the plate normal. This can be inferred from the photograph in Fig. 5 and seen more clearly in Fig. 7, which illustrates a view along the ScanHead® axis. Suitable angles for the stimulation and collection sections were determined by balancing the optical behavior of the incident beam (flare, scatter), the signal output at the photodetector, and MTF.

## 2.2 Light-collection optics

Once stimulated, the IP emits (blue) light in proportion to the stored x-ray exposure. This light emission is roughly Lambertian, that is, it is emitted in all directions. The purpose of the light-collection optics is to collect and direct as much of this emitted light as possible onto the active surface of the photodetector. In flying-spot scanners, this is usually done with acrylic light pipes, fiber optics, mirrors or integrating cavities. As in a flying-spot scanner, maximum light collection efficiency is realized when the entrance surface of the collection optics is as close to the IP surface as possible (i.e., the numerical aperture is large). For the same reason, the exit surface must be (at least, optically) close to the photodetector. The size constraints of ScanHead® add additional complexity, as already indicated. This produces a configuration for the reflective scanning geometry as shown in Fig. 7.

One solution for light collection involves an array of tiny gradient-index (GRIN) lenses, e.g., SELFOC® microlenses<sup>13</sup>, positioned along the illuminated line. These tiny (~1 mm diameter) lenses act like focused fiber optics that collect the light emitted from the scan line and image it, at unit magnification, onto the individual elements of the photodetector. In order to capture enough of the light emitted from the stimulated line, this SELFOC® array must contain multiple rows of hexagonally close-packed GRIN lenses across the entire 43 cm width of the scan line. In this configuration, the collection efficiency of the SELFOC® array is about 10%. This relatively low value is partly due to the empty interspaces between the lenses. In addition, GRIN optics are fairly dispersive for the main wavelengths (~ 400-450 nm) emitted by the storage phosphors of interest.

Another solution is the use of linear microlens arrays. The light transmission characteristics of such arrays can be significantly better (~2x) than GRIN optics, albeit at significantly higher cost. This is the solution chosen for ScanHead®. Fig.8 shows a small section of one possible microlens configuration for the "Collection Optics" block from Fig. 7. The array of microlenses is clearly visible between the two axial cylindrical lenses. The light collection efficiency in the direction perpendicular to the scan line is determined essentially by the height of the axial cylindrical lenses. Along the scan line, the efficiency is determined primarily by the width of each microlens. With appropriate angulation and spacing of the collection optics relative to the IP, the collection efficiency of the microlens array in ScanHead® is greater than 20%, comparable to that of light collection systems in flying-spot scanners.

## 2.3 Optical Filter

The purpose of the optical filter is to block all light that is not stimulated emission from the IP. This is critical since the ratio between the stimulating (red) light intensity and the emitted (blue) light intensity is greater than  $10^8$ . This means that the optical filter must have high transmission for the emitted wavelengths (400-450 nm, depending on the storage phosphor) and extremely low transmission for all other wavelengths to which the photodetector is sensitive.

The design of the optical filters for ScanHead® is more difficult than for all current flying-spot scanners. For one, current flying-spot scanners are designed to scan only one storage phosphor material (either BaFBr:Eu<sup>2+</sup> or RbBr:Tl). ScanHead®, on the other hand, will be configurable for two storage phosphor materials, BaFBr:Eu<sup>2+</sup> and the new CsBr:Eu<sup>2+</sup> needle-based phosphors, that have different emission characteristics<sup>3</sup>. Secondly, IR radiation plays a more significant role in ScanHead® (for example, the laser diode stimulation source produces IR radiation in addition to the visible). IR radiation is not a problem in laser-diode-based flying-spot scanners using photomultiplier tubes (PMTs) as photodetectors, since their response decreases rapidly above 650 nm. However, the CCD photodetectors used in ScanHead® are still quite sensitive in the IR, so the filter solution must block these wavelengths in addition to the stimulating light.

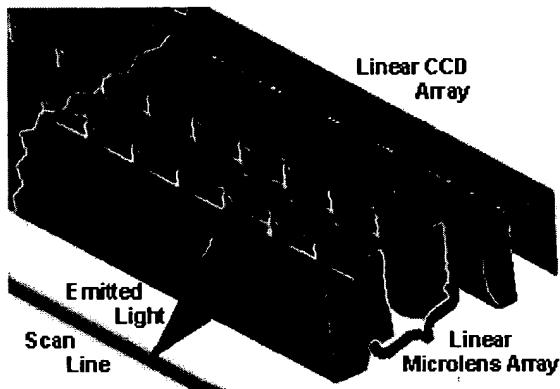


Figure 8. Linear microlens array used in the collection optics of ScanHead®. Behind the collection optics (optical filter not shown) is the CCD photodetector.

## 2.4 Photodetector and on-board electronics

The photodetector in ScanHead® (see Figs. 7 and 8) is an array of six custom-designed, low-noise, linear CCD sensors with on-board processing electronics and built-in temperature sensors. The six CCDs are butted end-to-end in a special housing to provide a practically continuous line of detector elements along the illuminated scan line (dead zones no larger than 70 $\mu\text{m}$ ). Each CCD contains 1464 asymmetric, active photo-elements measuring 50  $\mu\text{m} \times 400 \mu\text{m}$ , with the small dimension oriented along the illuminated scan line (the "fast-scan" direction in a flying-spot scanner). The asymmetry of the active elements allows the system to maintain high resolution along the scan line while still collecting enough emitted light per pixel to keep the system gain and SNR at reasonable levels (the resolution in the slow-scan direction is determined to first order by the width of the illuminated line and by light scattering within the phosphor layer). The signals from the 8784 individual elements (pixels) along the line can be combined in several different ways, leading to in-line sampling resolutions of 50, 100 or 150  $\mu\text{m}$  per pixel.

With appropriate driving circuitry and operating points, the dynamic range of the ScanHead® CCDs is over 10<sup>4</sup>, comparable to the PMTs used in flying-spot scanners, and well matched to the dynamic range of the storage phosphor screens used in diagnostic imaging. Thus, no special pre-scanning or gain setting is required to determine the actual exposure range on an IP. The quantum efficiency of the photodetectors for the emission wavelengths of interest is greater than 60%, which is better than that of PMTs in the same wavelength range.

The on-board electronics in ScanHead® take the analog video signal emerging from the CCDs, sample it and quantize it to 14 bits. Special effort was made in the board designs to keep electronic noise to a minimum. For example, since the internal noise of CCDs is normally greater than that of PMTs, the A/D converter uses correlated double sampling<sup>14</sup> to improve the read-out SNR. In this scheme, both the total signal value and the dark current (offset) are measured for each pixel, and then subtracted to produce the signal value that is quantized.

## 2.5 Other Scanner Components


While the ScanHead® is novel and interesting technology in and of itself, it is useless without the appropriate hardware and software that enable construction of a complete CR system. For example, the ScanHead® module must be moved at constant velocity along the IP (slow-scan transport). It must also maintain a constant height above the IP. Even small fluctuations in velocity or height can lead to visible banding in the image. The IP must be positioned and affixed to the scanning stage during ScanHead®'s movement. Furthermore, the scanned IP must be erased before re-use. This requires a new, integrated erase unit in a size consistent with ScanHead®'s (erase modules in current flying-spot systems are too large for ScanHead®-based systems). Solutions for some of these issues will be covered in more detail in forthcoming publications.

## 3. IMAGE QUALITY EVALUATION

One important measure of the imaging performance of a new capture system is its objective image quality relative to other systems designed for the same purpose. A suitable metric for such comparisons is Detective Quantum Efficiency. Here we must distinguish between the DQE results achieved with conventional powder IPs, where the phosphor layer consists of small phosphor particles embedded in a binder, and those achieved with needle IPs, where the phosphor layer is grown in tiny, vertical needles on a substrate, without the need for a binding material. The image quality results for the two are quite different<sup>3</sup>.


Shown in Fig. 9a are the DQE results for three different CR systems: a state-of-the-art flying-spot scanner reading a BaFBr<sub>x</sub>:Eu<sup>2+</sup> powder IP (Agfa ADC Compact/MD-30 screen), a prototype ScanHead®-based reader scanning the same (MD-30) powder IP, and the same reader scanning a prototype CsBr:Eu<sup>2+</sup> needle IP. All data were acquired at 70 kV<sub>p</sub> with a 0.5 mm Cu filter at a dose of 1  $\mu\text{Gy}$ . The Scan Head® collection optics for both screens were fitted with an early prototype linear microlens array and, to illustrate the impact of this component on image quality, both the in-line ("fast-scan") and slow-scan DQE results are shown for the ScanHead®/Needle IP configuration (see below). The experimental procedures used to measure the MTF and NPS, and the methodology to calculate DQE in the CR systems have been described in more detail elsewhere<sup>3</sup>.

As can be seen in the figure, the combination of ScanHead® (with a microlens) and a powder IP produces DQE results that are comparable to, or somewhat better than today's flying-spot systems scanning powder IPs. The DQE results for

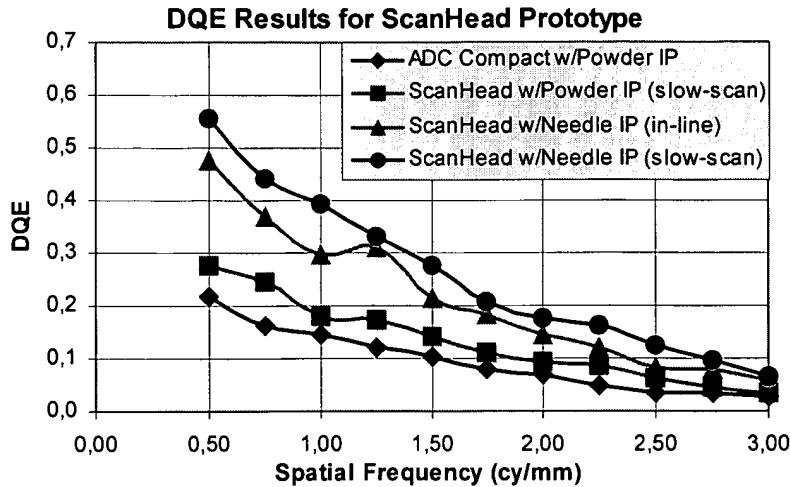

**Figure 9a.** DQE results for three imaging systems: CR using a flying-spot scanner (ADC Compact) and a conventional (MD-30) powder IP, CR using a ScanHead® prototype and a conventional (MD-30) powder IP (only slow-scan DQE shown), CR using a ScanHead® prototype with a prototype CsBr:Eu<sup>2+</sup> needle IP (in-line and slow-scan DQEs are shown). The difference between in-line and slow-scan DQE for this ScanHead® system is due to the linear microlens used for these measurements (see text). All data were acquired at a dose of 1  $\mu\text{Gy}$  at 70 kV<sub>p</sub> with a 0.5 mm Cu filter.

ScanHead® and Needle IP are significantly better (2-3x) than those produced with powder IPs. This improvement is due largely to the superior imaging properties of the new needle screens<sup>3</sup> but also to the imaging performance of ScanHead®. Note that the slow-scan and in-line DQE for this ScanHead® module are different. The linear microlens array used in this prototype scanner had suboptimum MTF characteristics along the in-line direction, causing a decrease in DQE. More recent versions of this array have shown better optical properties, and we expect the in-line results in future ScanHead® units to be equivalent to those in the slow-scan direction.

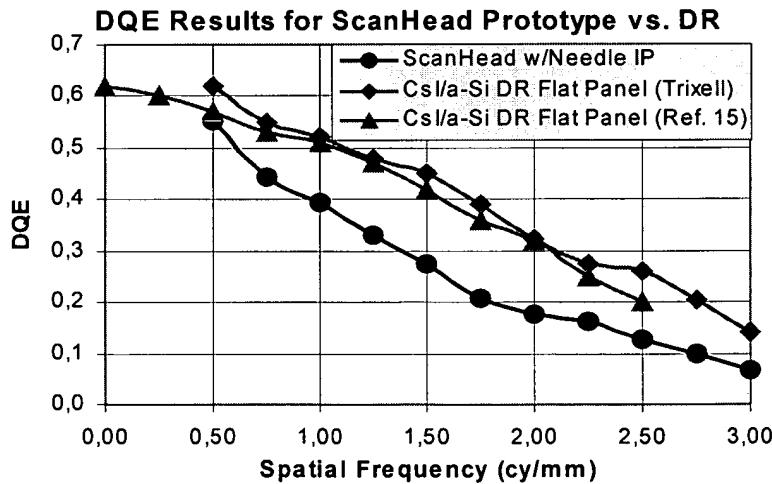


**Figure 9b.** DQE results (at 1  $\mu\text{Gy}$ ) for a ScanHead® prototype with a prototype needle IP, along with two flat-panel DR systems based on CsI and an amorphous Si photodetector/converter<sup>15</sup>. The DQE data for the Trixell detector (◆) were acquired in-house using the same methodology as in the CR measurements shown in Fig. 9a.

It is also interesting to compare the ScanHead®/Needle IP DQE results with those of DR systems. Unfortunately, DQE comparisons between systems can be difficult and dangerous to interpret due to differences in exposure, measurement and data analysis techniques. With this caveat, Fig. 9b shows a comparison of the ScanHead®/Needle IP results from Fig. 9a (slow-scan) and literature values<sup>15</sup> on one flat-panel DR system based on a CsI:Tl x-ray detector and an amorphous Si photodetector array. Another DR flat panel (from Trixell) was measured on an in-house system with the same methodology used for the three CR systems shown in Fig. 9a. These results are also shown.


Although there is currently a gap (30-50%) between the ScanHead® and DR DQE curves in Fig. 9b, particularly at mid-frequencies, it is important to note that the CR data were created on a laboratory version of the device using early prototype needle screens. We expect the ScanHead® DQE values to increase as the needle IPs improve (increase in gain and sharpness) and as the optical filters improve (decrease in electronic noise). Despite the current gap, the visual image quality of the ScanHead®/Needle images is comparable to that of DR systems. A recent contrast-detail study comparing a conventional flying-spot scanner using powder IPs, a ScanHead® scanner using Needle IPs, and several flat-panel DR systems suggests that this objective image quality parity carries over into observer performance for detection tasks with low-contrast targets in noisy backgrounds<sup>12</sup>.

#### 4. APPLICATIONS AND DEVICES

A compact scanning technology like ScanHead®, with its high level of image quality, can be used in practically any diagnostic imaging application and in multiple device configurations. The simplest configuration is a single-slot, cassette-based scanner for distributed or low-volume imaging applications (e.g., in the exam room). One possible design for such a scanner is illustrated in Fig. 10. Note that the scanning engine (ScanHead®), its slow-scan transport system, an erase unit, the IP handling/positioning unit and platform, and all associated electromechanical and electronic systems fit inside a floor-standing scanner "tower." The scanner can be adapted to read conventional powder IPs or needle IPs, assuming the same cassette design. With appropriate cassette buffering, centralized, high-volume configurations can also be realized with the same basic components.

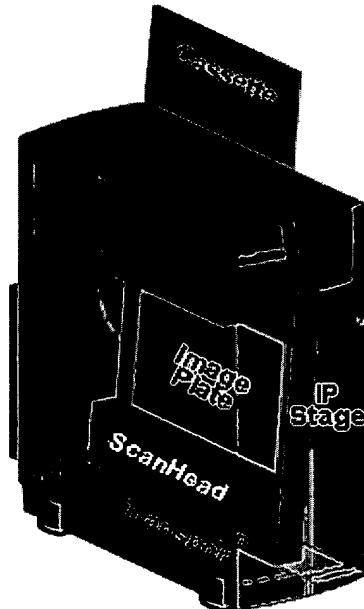


**Figure 10.** A view through the semi-transparent case of a single-slot CR tower reader based on ScanHead® technology. The tower contains the scanning engine, the scanning stage, slow-scan transport hardware, an erase unit, IP handling and positioning subsystems, and all associated electronics.

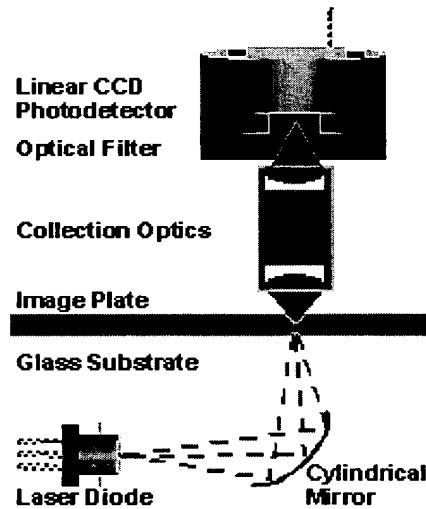

**Figure 11.** Possible transmission-mode scanning geometry using a cylindrical mirror. The stimulation section of ScanHead® is below, the collection section above the IP.

ScanHead® technology also enables some devices that were previously impossible with flying-spot CR technologies. For example, flat-panel, cassetteless CR systems can now be realized. For these devices, a transmissive ScanHead® configuration is more useful than the reflective one discussed so far. In this configuration, the stimulation and collection sections are on opposite sides of the IP. Fig. 11 shows a possible beam geometry of such a system and Fig. 12 shows a prototype of a ScanHead®-based flat-panel CR system that might be used for portable radiography or as the capture

module inside a Bucky table. For flat-panel systems, the IP must be optically accessible from both sides, so the substrate must be transparent, e.g., glass. Since the needle-based storage phosphor, CsBr:Eu<sup>2+</sup> can be deposited very effectively on glass substrates, and produces excellent image quality, it is the ideal candidate for the CR flat-panel. The radiation sensitivity of the CCDs can be dealt with by parking the ScanHead® module in a shielded area during the x-ray exposure. One can also imagine other possible ScanHead® device configurations suitable for a dedicated chest unit, mammography systems, etc.

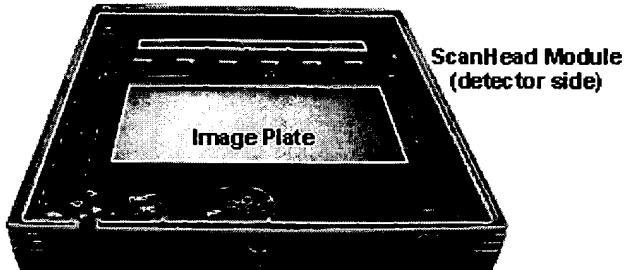

**Figure 12.** View inside a CR flat panel containing the ScanHead® module. The stimulation section is underneath the image plate.

## 5. CONCLUSIONS

We have described a new type of CR scanning engine using a **line-at-a-time** scanning principle, instead of the point-at-a-time technique used by all current CR systems. This produces a number of significant advantages over current flying-spot scanners. The size of scanners using this ScanHead® technology can be considerably smaller than current systems due to the fact that the stimulation and collection subsystems are contained in one compact module. Despite the size reduction, scanner throughput can still be quite high since many pixels are captured simultaneously.

The image quality of a ScanHead® system can be equivalent to current CR systems when using conventional powder IPs, and can be significantly better than current CR systems when using needle IPs. The high level of image quality achieved with needle IPs is largely due to the properties of the needle IP itself, but also due to the imaging characteristics of ScanHead® (longer pixel integration time, higher read-out depth, good MTF, low-noise electronics). With appropriate collection optics, the image quality with needle IPs can be comparable to that of recently introduced flat-panel DR systems based on CsI:Tl and amorphous Si photodetectors. While image quality is certainly not the only criterion for the clinical applicability or acceptability of a new technology, it is not an insignificant one. The variety of device configurations possible with ScanHead® technology, from centralized, high-volume systems to distributed, single-plate readers to cassette-less flat panels, and its significantly lower cost compared to current DR systems, suggest that future CR systems will continue to play an important role in the acquisition of diagnostic images.

## ACKNOWLEDGEMENTS

We wish to thank Magath Niang for performing the DQE measurements. We also wish to thank Dr. Georg Reiser for helpful conversations on ScanHead® optics. Thanks also to Christian Bommer for help with the 3-d microlens graphic.

## REFERENCES

1. R. Schaetzing, B. R. Whiting, A. R. Lubinsky, J. F. Owen, "Digital Radiography Using Storage Phosphors," *Digital Imaging in Diagnostic Radiology*, eds. J. D. Newell, C. A. Kelsey, pp. 107-138, Churchill Livingstone, New York, 1990.
2. M. Sonoda, M. Takano, J. Miyahara, H. Kato, "Computed Radiography Utilizing Scanning Laser Stimulated Luminescence," *Radiology* **148**:3, pp. 833-838, 1983.
3. P. Leblans, L. Struye, "New Needle-crystalline CR Detector," *Proc. SPIE* **4320**, pp. 59-67, 2001.
4. H. von Seggern, T. Voigt, K. Schwarzmichel, "Physical Mechanisms of Photostimulation in Medical X-Ray Storage Phosphors," *Siemens Forsch.-u. Entwickl.-Ber.* **Bd 17**:3, p.125, 1988.

5. H. von Seggern, "Photostimulierbare Speicherleuchtstoffe für Röntgenstrahlung," *Phys. Bl.* **48**:9, pp. 719-723, 1992.
6. P. Leblans, private communication.
7. M. I Sezan, A. M. Tekalp, "Survey of Recent Developments in Digital Image Restoration," *Opt. Eng.* **29**(5), pp. 393-404, 1990.
8. J. C. Dainty, R. Shaw, *Image Science*, Academic Press, London, 1976.
9. H. Kato, "Photostimulable Phosphor Radiography Design Considerations," Proc. AAPM Summer School, pp. 860-898, 1991.
10. J. T. Dobbins III, D. L. Ergun, L. Rutz, D. A. Hinshaw, H. Blume, D. C. Clark, "DQE(f) of Four Generations of Computed Radiography Acquisition Devices," *Med. Phys.* **22**(10), pp. 1581-1593, 1995.
11. S. Arakawa, W. Itoh, K. Kohda, T. Suzuki, "Novel Computed Radiography System with Improved Image Quality by Detection of Emissions from Both Sides of an Imaging Plate," *Proc. SPIE* **3659**, pp. 572-581, 1999.
12. R. Schaetzing, R. Brandl, K.-J. Pfeifer, P. Leblans, "Improving the DQE and Low-Contrast Detectability of CR Systems Through the Use of Structured (Needle) Storage-Phosphor Screens," *Radiology* **221**(P), 463, 2001.
13. *SEFOC* is a registered trademark of Nippon Sheet Glass Co., Ltd., Japan.
14. Specifications Datasheet for Analog Devices AD9814 CCD/CIS Signal Processor, Rev. 0, Analog Devices Inc., Norwood, MA, 1999.
15. K. S. Kump, "Fast Imaging of a 41cm Amorphous Silicon Flat-Panel Detector for Radiographic Applications," *Proc. SPIE* **4320**, pp. 87-93, 2001.


**POERSCH RECIPROCATING BUCKYS, MODELS RRJ & RRK  
SERVICING AND REPLACEMENT PARTS**

# 115 VOLTS A.C.


TERMINAL MODULE  
PIN 8284R19




REPLACEMENT FUSES  
QUAN. ONE FOD (2)PFS  
TIME DELAY TYPE FST  
5 x 20 mm, 500 mA, 250 Volts  
PIN 8284R19-01

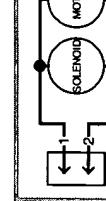


DRIVE MODULE  
PIN 8284R49

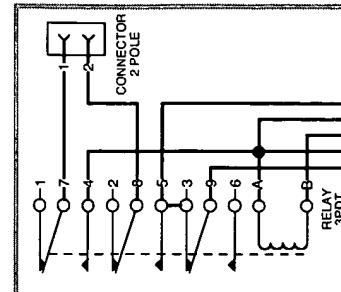

RELAY  
PIN 8284R29


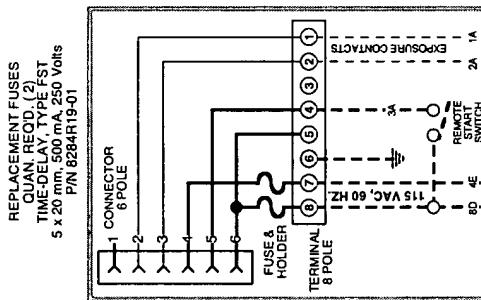


ALTERNATE 5 WIRE HOOK-UP  
WITH EXTERNAL FUSING

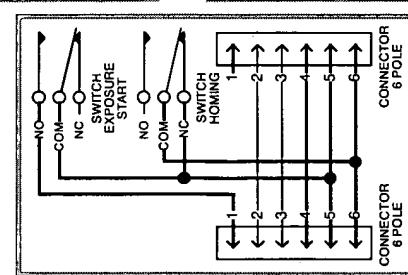


SWITCH MODULE  
PIN 8284R39

RELAY  
PIN 8284R29

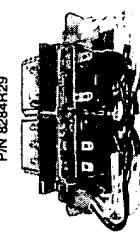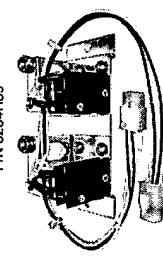


## SCHEMATIC WIRING DIAGRAM


LEGEND	115 VAC Wiring, By Manufacturer
- - -	115 VAC Wiring, By Others
—	Low-Voltage Wiring, By Manufacturer
— - -	Low-Voltage Wiring, By Others

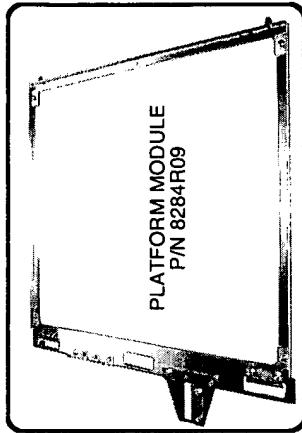


PLATFORM MODULE  
PIN 8284R09

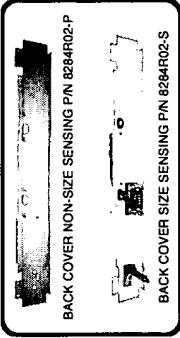

BACK COVER NON-SIZE SENSING PIN 8284R02-S

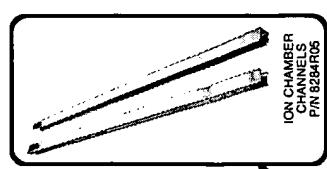


ION CHAMBERS  
PIN 8284R05

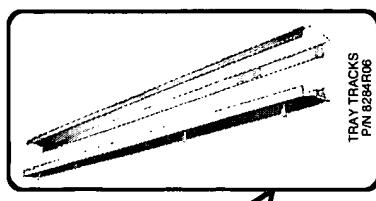

TRAY TRACKS  
PIN 8284R06



TRAY STOP  
PIN 8284R04




DRIVE MODULE  
PIN 8284R19

When ordering replacement parts provide MODEL AND SERIAL NUMBER  
For service and parts contact:

Poersch Metal Mfg. Co.  
4027 W. Kinzie St.  
Chicago, IL 60624 USA  
Phone: 773.722.0890, Fax: 773.722.4122  
Web: www.poerschmetal.com

## INSTRUCTIONS HOW TO MANUALLY MOVE RECIPROCATING PLATFORM

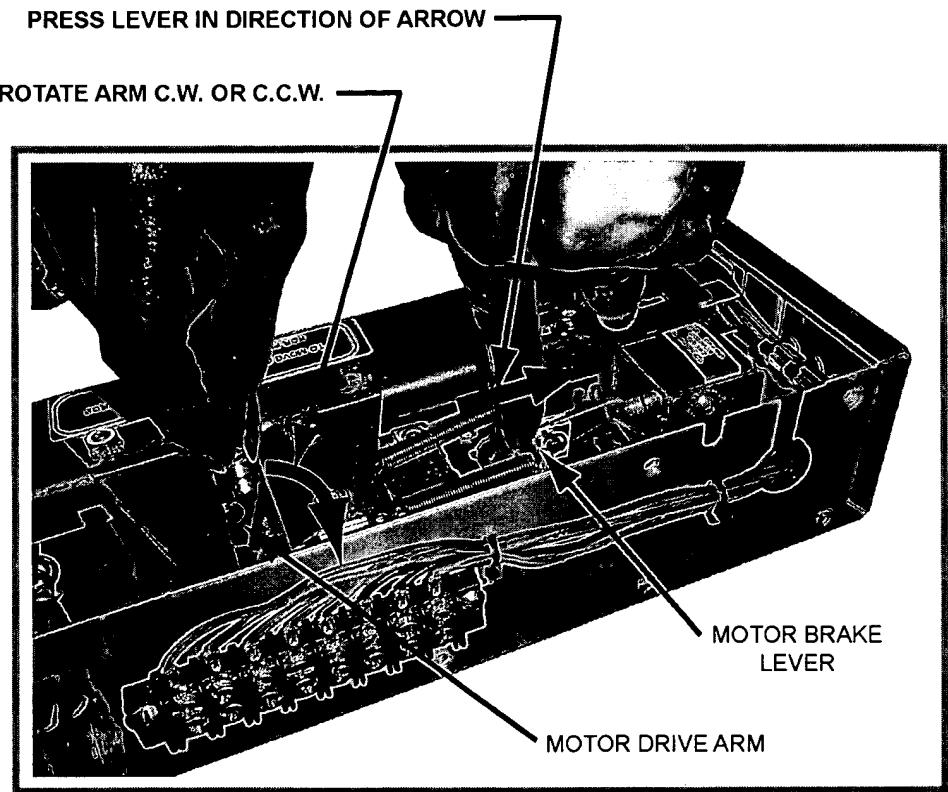


**STEP 1** DISCONNECT POWER TO BUCKY.

**STEP 2** Remove bucky motor compartment top cover.


**STEP 3** Release motor brake by pressing lever back with finger, in direction of arrow, away from springs. (see photo above)

**STEP 4** Rotate motor drive arm in either direction. (see photo above)  
Note, bucky will automatically reset itself when power is turn on.

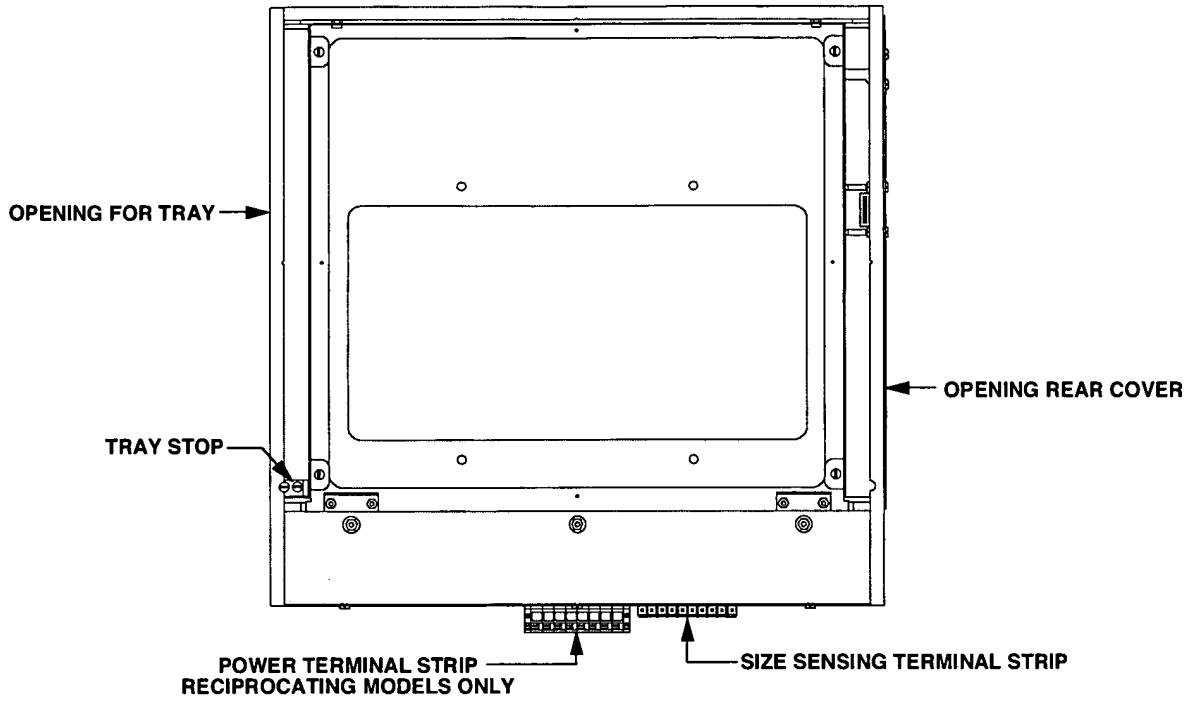

### BUCKY - LEFT HAND LOADING

TURN PAGE FOR INSTRUCTIONS  
HOW TO SWITCH LOADING SIDE FROM LEFT TO RIGHT

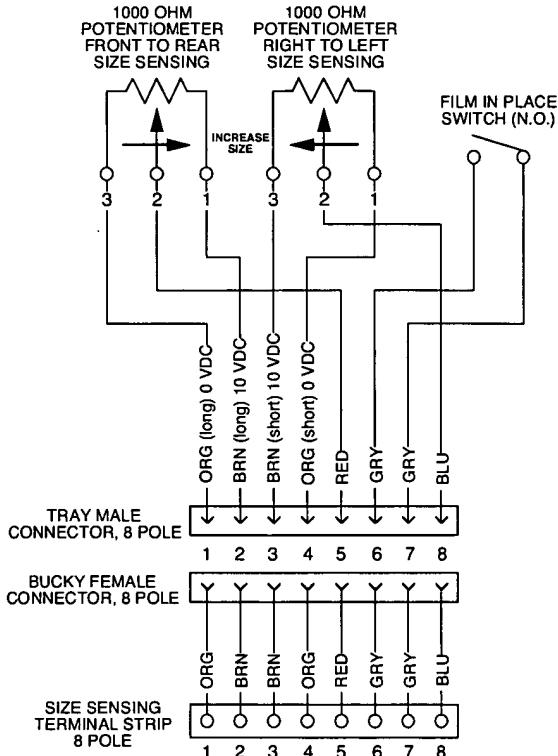


### OUTPUT VOLTAGE


With  $10 \pm 0.01$  VDC applied across potentiometers (pins 1 and 2 for front to rear sizing; pins 3 and 4 for right to left sizing), the voltage measured between terminals 1 and 5 (for front to rear sizing) and 4 and 8 (for right to left sizing) for various cassette sizes listed below.

CASSETTE SIZE (OUTSIDE DIM.)	STATED FILM SIZE	OUTPUT VOLTAGE $\pm 0.2$	OHMS VALUES
154.5 mm	5 in.	0.60	940
157.5 mm	13 cm	0.69	931
205.3 mm	7 in.	2.06	794
207.5 mm	18 cm	2.12	788
230.7 mm	8 in.	2.80	720
268.8 mm	24 cm	3.90	610
281.5 mm	10 in.	4.26	574
306.9 mm	11 in.	5.00	500
327.5 mm	30 cm	5.60	440
332.3 mm	12 in.	5.70	430
383.1 mm	35 cm	7.02	298
459.3 mm	43 cm	9.40	60


WIRING SCHEMATIC FOR SIZE SENSING  
MODEL RRJ BUCKYS, MODEL RSJ GRID CABINETS  
AND MODEL QJ TRAYS

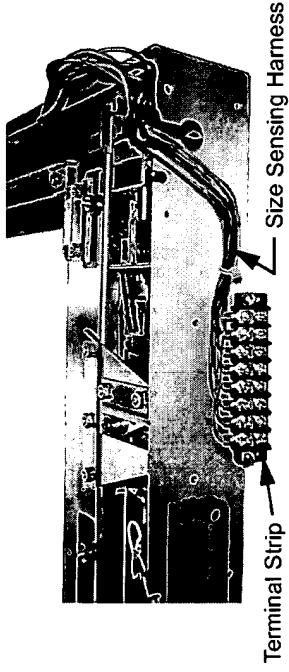


**STEP 1**  
Remove bucky motor compartment top cover and unscrew size sensing terminal strip from outside of bucky (Note: do NOT disconnect wires from terminal strip). Cut two wire ties binding size sensing harness inside motor compartment. Remove grommet and lift out harness from hole on side of bucky.

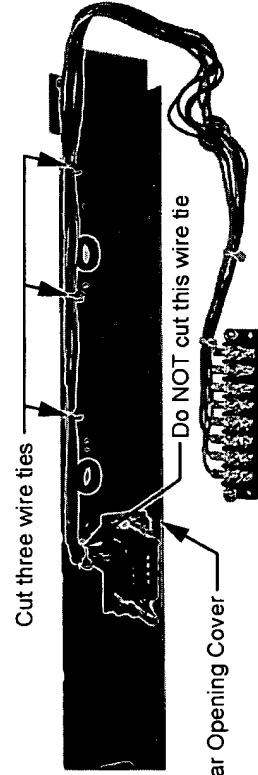


**STEP 2**  
Unscrew four screws fastening tray opening, rear cover to bucky, remove rear cover from bucky. Cut three wire ties binding harness along top flange of rear cover. Do NOT cut wire tie above the blue female connector assembly.

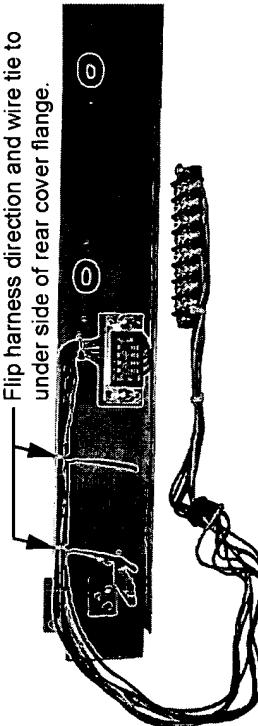


**STEP 3**  
Flip size sensing harness 180° to other direction and wire tie securely along under side of rear cover top flange.

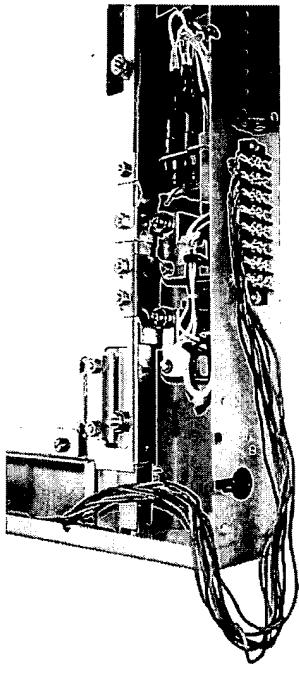

**STEP 4**  
Unscrew "tray stop" from bottom of tray side of bucky; refasten stop at opposite diagonal corner. Reinstall rear cover into opposite opening of bucky and fasten size sensing terminal strip to outside of bucky.

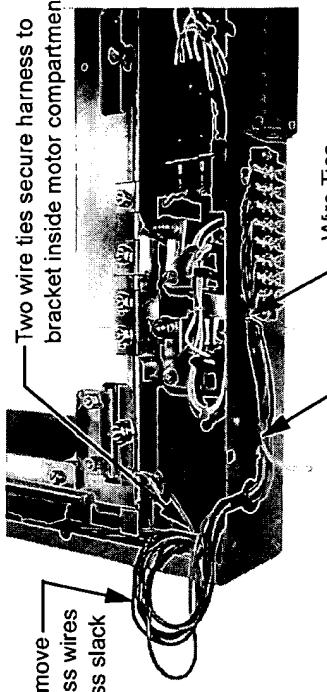

**STEP 5**  
Remove grommet and dress harness wires thru hole on side of bucky, reinsert grommet. Secure harness with two wire ties inside motor compartment and two wire ties outside bucky. BEFORE pulling wire ties tight remove slack from each wire of harness as shown above.

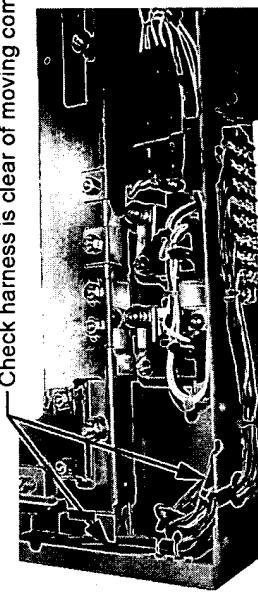


**STEP 6**  
Tightly bundle harness excess wire and secure inside motor compartment away from all moving components. Replace and fasten top cover securely.

## INSTRUCTIONS HOW TO SWITCH BUCKY LOADING SIDE FROM LEFT TO RIGHT

FOLLOW SIMILAR PROCEDURE WHEN SWITCHING FROM RIGHT TO LEFT

**POERSCH** Metal Manufacturing Co.  
4027 WEST KINZIE STREET, CHICAGO, IL 60624, USA  
PHONE: 773.722.0890; FAX: 773.722.4122  
[www.poerschmetal.com](http://www.poerschmetal.com)

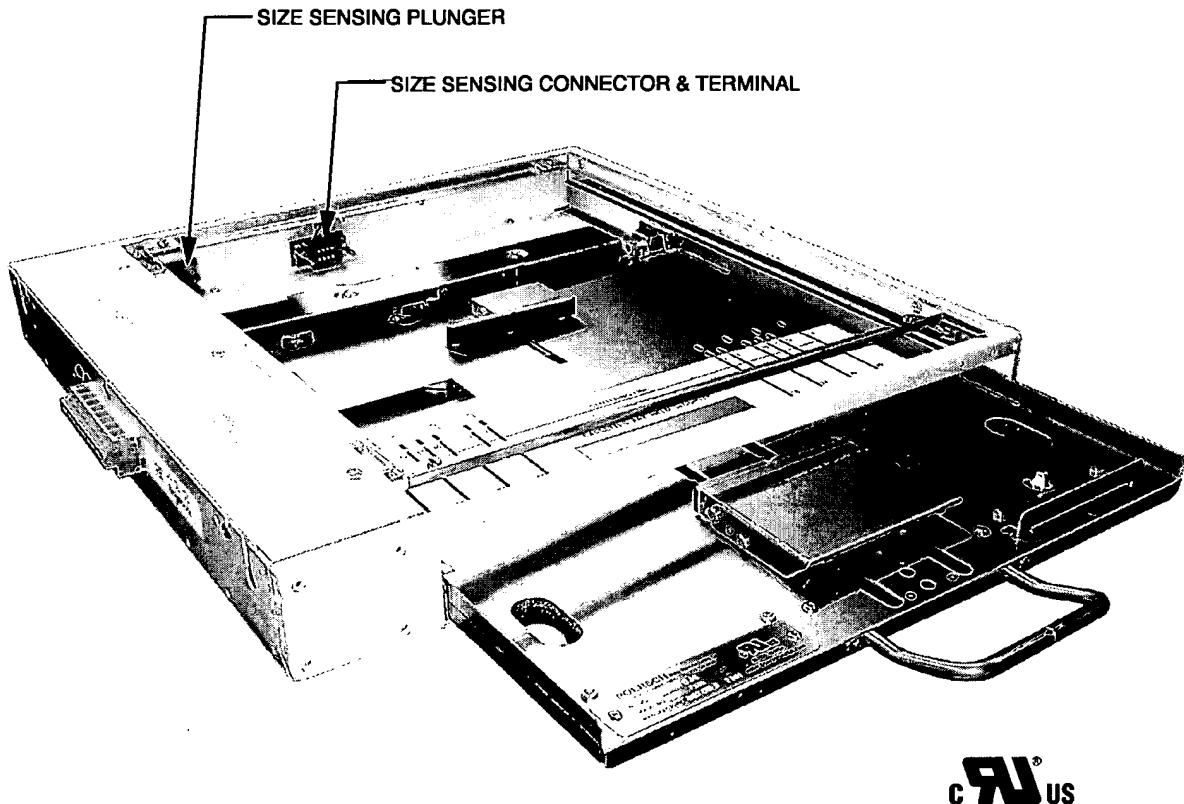

### RECIPROCATING BUCKY MODEL RRJ


Shown with automatic size sensing cassette tray partially inserted  
(cassette tray not included)

DESCRIPTION	PART NO.
AUTOMATIC SIZE SENSING INSTALLATIONS	RRJ
MANUAL SIZE SENSING INSTALLATIONS	RRK

The Poersch reciprocating bucky is designed for use in Radiographic Tables and Vertical Bucky Stands. Available in two models; model RRJ is compatible with automatic size sensing cassette trays and model RRK is compatible with manual non-size sensing cassette trays. Both models are identical in construction and size. The plunger, connector and terminal required for automatic size sensing are not included on the model RRK.

The cassette loading, is easily switched from right hand to left hand by transferring the rear cover and tray stop to opposite ends. Under normal conditions no maintenance is required. All drive components are rated for continuous duty and life tested for greater than one million exposures.

**GENERAL SPECIFICATIONS**  
**FINISH:** EXCEPT WHERE SPECIFIED, SPRAYLAT  
 POWDER COATING, PAPYRUS WHITE  
 HYBRID TEXTURE, MATTE  
**TOLERANCES:** .000 = ± .007", .000 = ± .004"  
**WEIGHTS:** NET 22 LBS., SHIPPING 28 LBS.  
 VOLTAGE: 120 VOLTS, 60 HZ.  
**AMPERAGE:** 0.5 AMPS

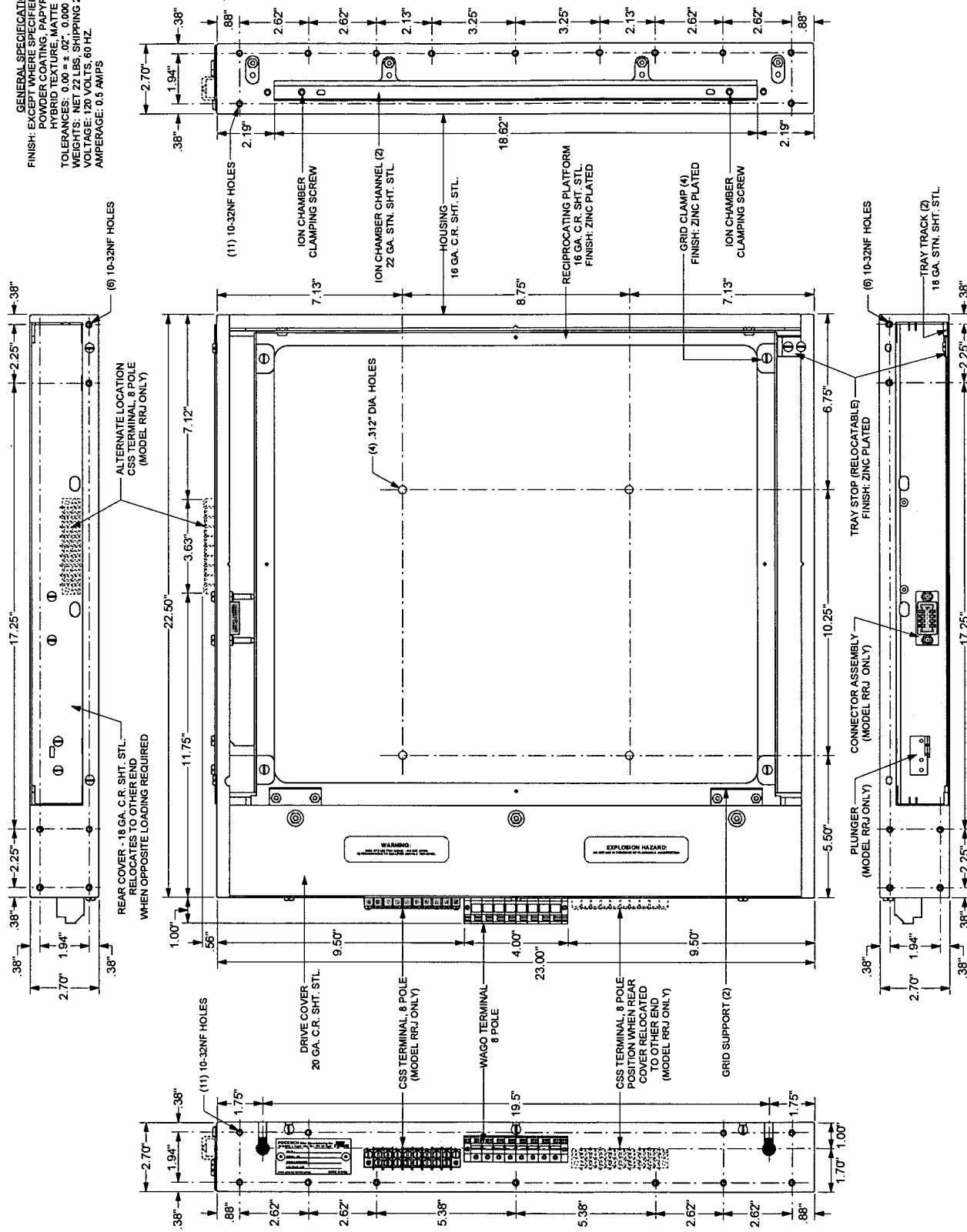

## BUCKY (RECIPROCATING) MODEL RRJ AND RRK (RIGHT HAND LOADING ILLUSTRATED)